



INVESTOR IN PEOPLE



The Patent Office
Concept House
Cardiff Road
Newport
South Wales
NP10 8QQ

The undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in the certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

Signed

Dated

- 7 DEC 2000

Patents Act 1977
(Rule 16)

Request for grant of a patent

(See the notes on the back of this form. You can also get an explanatory leaflet from the Patent Office to help you fill in this form)



23NOV99 E493801-4 D00027
P01/7700 0.00-9927590.1

The Patent Office
Cardiff Road
Newport
Gwent NP9 1RH

1. Your reference	85.71394		
2. Patent application number (The Patent Office will fill in this part)	9927590.1		22 NOV 1999
3. Full name, address and postcode of the or of each applicant (<i>underline all surnames</i>)	GENPAK Limited Science Park Square Falmer, Brighton BN1 9SB Sussex, U.K.		
Patents ADP number (<i>if you know it</i>) 7704331001			
If the applicant is a corporate body, give country/state of incorporation			
4. Title of the invention	Microarray slides		
5. Name of your agent (<i>if you have one</i>)	Frank B. Dehn & Co.		
"Address for service" in the United Kingdom to which all correspondence should be sent (including the postcode)	179 Queen Victoria Street London EC4V 4EL		
Patents ADP number (<i>if you know it</i>)	166001 ✓		
6. If you are declaring priority from one or more earlier patent applications, give the country and the date of filing of the or of each of these earlier applications and (<i>if you know it</i>) the or each application number	Country	Priority application number (<i>if you know it</i>)	Date of filing (day / month / year)
	—	—	—
7. If this application is divided or otherwise derived from an earlier UK application, give the number and the filing date of the earlier application	Number of earlier application	Date of filing (day / month / year)	
8. Is a statement of inventorship and of right to grant of a patent required in support of this request? (<i>Answer 'Yes' if:</i> a) any applicant named in part 3 is not an inventor, or b) there is an inventor who is not named as an applicant, or c) any named applicant is a corporate body. See note (d))	Yes		

Patents Form 1/77

9. Enter the number of sheets for any of the following items you are filing with this form. Do not count copies of the same document

Continuation sheets of this form

-

Description

3

Claim(s)

4

Abstract

1

Drawing(s)

1 + 1

10. If you are also filing any of the following, state how many against each item.

Priority documents

-

Translations of priority documents

-

Statement of inventorship and right to grant of a patent (*Patents Form 7/77*)

-

Request for preliminary examination and search (*Patents Form 9/77*)

1

Request for substantive examination (*Patents Form 10/77*)

1

Any other documents
(*please specify*)

-

11. I/We request the grant of a patent on the basis of this application.

Signature

Date 19 November 1999

12. Name and daytime telephone number of person to contact in the United Kingdom

P.M. Jeffrey
01273 244200

Warning

After an application for a patent has been filed, the Comptroller of the Patent Office will consider whether publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore, if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a patent abroad without first getting written permission from the Patent Office unless an application has been filed at least 6 weeks beforehand in the United Kingdom for a patent for the same invention and either no direction prohibiting publication or communication has been given, or any such direction has been revoked.

Notes

- If you need help to fill in this form or you have any questions, please contact the Patent Office on 0645 500505.
- Write your answers in capital letters using black ink or you may type them.
- If there is not enough space for all the relevant details on any part of this form, please continue on a separate sheet of paper and write "see continuation sheet" in the relevant part(s) of the form. Any continuation sheet should be attached to this form.
- If you have answered 'Yes', Patents Form 7/77 will need to be filed.
- Once you have filled in the form you must remember to sign and date it.
- For details of the fee and ways to pay please contact the Patent Office.

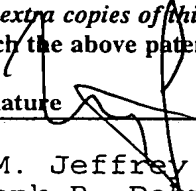


The
Patent
Office

7/77

Statement of inventorship and of
right to grant of a patent

The Patent Office
Cardiff Road
Newport
Gwent NP9 1RH

1.	Your reference	85.71394
2.	Patent application number (if you know it)	9927590.1
3.	Full name of the or of each applicant	GENPAK Limited
4.	Title of the invention	Adhesive Label with grid for microscope slide
5.	State how the applicant(s) derived the right from the inventor(s) to be granted a patent	By virtue of contracts of employment
6.	How many, if any, additional Patents Forms 7/77 are attached to this form? (see note (c))	-
7.	<p>I/We believe that the person(s) named over the page (and on any extra copies of this form) is/are the inventor(s) of the invention which the above patent application relates to.</p> <p>Signature  Date 21 November 2000</p>	
8.	Name and daytime telephone number of person to contact in the United Kingdom	<p>P.M. Jeffrey Frank B. Dahn & Co. 01273 244200</p>

Notes

- a) If you need help to fill in this form or you have any questions, please contact the Patent Office on 0645 500505.
- b) Write your answers in capital letters using black ink or you may type them.
- c) If there any more than three inventors, please write the names and addresses of the other inventors on the back of another Patents Form 7/77 and attach it to this form.
- d) When an application does not declare any priority, or declares priority from an earlier UK application, you must provide enough copies of this form so that the Patent Office can send one to each inventor who is not an applicant.
- e) Once you have filled in the form you must remember to sign and date it.

Enter the full names, addresses and postcodes of the inventors in the boxes and underline the surnames

Julian BURKE
1 Apsley Cottage
Alciston
Polegate
East Sussex
BN26 6UR
UK

Patents ADP number (if you know it): 7955424CC1

Nishith PATEL
7 D'Aubigny Road
Brighton
BN2 9SB
UK

Patents ADP number (if you know it): 8029243CC1

Reminder

Have you signed the form?

Patents ADP number (if you know it):

71394.348

MICROARRAY SLIDES

5 The present invention relates to microarray slides, particularly microarray slides used in manual microarraying.

Microarraying involves the transfer of e.g. biological DNA material from a source e.g. a microtitre plate to a target e.g. a glass microscope or microarray
10 slide. Microtitre plates holding 96, 384 or 1534 different DNA samples are known, and conventionally it is desired to accurately transfer the samples from the source onto the target in the form of small distinct and separate micro spots. Once the material has been
15 deposited in the form of micro spots the individual samples can then be analysed.

The accurate transfer of material from a source to a target can be achieved by either robotic or manual spotting. For robotic transfer the configuration of the
20 material on the source can be easily monitored and controlled by setting the robotic system to spot in pre-defined areas and in a pre-defined configuration.

However, for manual transfer of material, the
25 placement of the material from the source onto the target tends to be more difficult to monitor and control. The material being transferred is often clear and thus once it is dry it can be very difficult to know where a spot has already been placed and therefore
30 exactly where the next spot should be placed on the target.

It is therefore desired to overcome the problems associated with conventional manual spotting techniques.

According to the present invention there is
35 provided the combination of a microscope slide and an adhesive label as claimed in claim 1. Other aspects of the invention are recited in the additional independent claims.

According to a preferred embodiment a pre-gridded microscope slide is provided which has a pre-printed alpha-numeric grid attached to the underside of the slide. The alpha-numeric grid appears through the glass
5 slide with the numbers and letters in the correct configuration i.e. write read up.

Various embodiments of the present invention will now be described, by way of example only, and with reference to the accompanying drawings in which:

10 Fig. 1 shows a microscope slide with a label attached thereto; and

Figs. 2(a) and 2(b) show two different labels.

With reference to Fig. 1, transfer of material onto a target microscope slide 1 is carried out from a source
15 plate to pre-defined positions on the top surface of the glass slide 1, using a pre-printed grid 2, preferably printed on a label 3, as a guide for placement of material to be spotted. Knowing which box or cell of the grid 2 has previously been spotted enables a user to
20 safely spot the next sample in the next available empty box or cell of the grid or array 2.

The format of the pre-printed grid 2 may vary, but a 8 x 12 array or grid corresponding with a 96 well microtitre plate format is preferred. According to such
25 an embodiment, 96 different DNA samples can be transferred from the source microtitre plate to the target microscope slide 1 with a one to one correspondence between the two of them.

Two different embodiments of label design are shown
30 in Figs. 2(a) and 2(b). The labels are shown enlarged. In the embodiment shown in Fig. 2(a) each array 2 is approximately 6 mm x 4 mm i.e. each cell is approximately 0.5 mm x 0.5 mm. In the embodiment shown in Fig. 2(b) each array 2 is approximately 12 mm x 8 mm
35 i.e. each cell is approximately 1 mm x 1 mm. The size and number of grids 2 on a single microscope slide 1 may vary depending on the amount of material that needs to be transferred. For the transfer of relatively large

sample amounts, the grids 2 can be made correspondingly larger so that the individual cells of the grid 2 can hold sufficient material.

5 After spotting, the pre-gridded slides 1 can be used in further analysis stages which may require the slides 1 to withstand extreme temperatures. The pre-gridded slides 1 can preferably withstand repeated freezing, including temperatures down to -20°C and repeated heating, up to 96°C over prolonged periods of
10 time. The slides 1 are also preferably resistant to corrosive chemicals and reagents.

The grids 2 have also been shown not to interfere with scanning of the glass slides 1 to detect fluorescent dyes, which is usually one of the final
15 stages of spotted material analysis.

Preferably, the printed grids can be removed from the glass microscope/arraying slides 1, after transfer of material has occurred.

71394.348

Claims

- 5 1. A microscope slide in combination with a self-adhesive label, said microscope slide having a front face upon which, in use, a sample to be analysed is deposited, and a rear face, wherein:
- 10 said label comprises indicia means so that when said label is attached, in use, to the rear face of said microscope slide, said indicia means is visible when viewed through said front face of said microscope slide.
- 15 2. The combination as claimed in claim 1, wherein said indicia means comprises an alpha-numeric grid.
- 20 3. A microscope slide having a front face and a rear face, and wherein said microscope slide comprises an alpha-numeric grid.
- 25 4. A microscope slide as claimed in claim 3, wherein said alpha-numeric grid is attachable to the rear of said microscope slide so that said grid is viewable through said front face.
- 30 5. An adhesive label having a plurality of alpha-numeric grids printed thereon.
- 35 6. A label as claimed in claim 5, wherein the characters forming said alpha-numeric grids are printed as a mirror image of conventional alpha-numeric.
7. A pre-gridded microscope slide having a grid printed on or attached to one side.
8. A microscope slide as claimed in claim 7, wherein said grid comprises alpha-numeric characters.

9. A microscope slide as claimed in claim 8, wherein at least some of said alpha-numeric characters are printed so that when viewed in a mirror numbers will appear as 0, 1, 2, 3 ... n and characters will appear as A, B, C ... Z.

10. Apparatus for assisting in transferring biological or similar material to a microscope slide, said apparatus comprising a grid which is attachable to a microscope slide.

11. A microscope slide having an integral grid for assisting in the transfer of specimens to a desired location on said slide.

12. A microscope slide having at least one grid attached thereto, said grid having 96 or more cells.

13. An adhesive label for assisting the spotting of a microscope slide, said label having on one side an adhesive layer and on an opposed side a grid together with mirror images of alpha-numeric characters.

14. An arraying slide in combination with an attachable pre-printed display, said display dividing in use said slide into a plurality of areas each having an area of 1 mm² or less.

15. A microscope slide having at least 96 indexed areas in combination with a microtitre plate having at least 96 sample cells.

16. A method of transferring DNA samples from a microtitre plate to a microscope slide, comprising the steps of:

providing a microtitre plate with a plurality of DNA samples arranged in an array of samples;
providing a microscope slide;

attaching a pre-printed label having at least one grid to said microscope slide; and

transferring DNA material from said microtitre plate to said microscope slide using the grid on said label as a guide.

17. A method of manually spotting a microscope slide, comprising the steps of:

attaching a gridded indicia means to said microscope slide; and

spotting said microscope slide using said indicia means as a guide.

18. A microscope slide having at least three 8 x 12 arrays within an area of 5 cm x 1.5 cm provided thereon.

19. A microscope slide having at least six 8 x 12 arrays within an area of 3 cm x 1.5 cm provided thereon.

20. A microscope slide as claimed in claim 18 or 19, wherein said arrays are printed on an adhesive label.

21. A microscope slide as claimed in claim 18, 19 or 20, wherein said arrays comprise alpha-numeric characters.

22. A microscope slide as claimed in claim 21, wherein said characters are printed as mirror images of conventional alpha-numeric characters.

23. A plurality of glass microscope slides, each said slide having indicator means which defines at least 50 areas or cells, each said area or cell being smaller than or equal to 1 mm².

24. A plurality of glass microscope slides as claimed in claim 23, wherein each said area or cell is smaller than or equal to 0.5 mm².

25. A plurality of glass microscope slides as claimed in claim 24, wherein each said area or cell is smaller than or equal to 0.25 mm^2 .

- 5 26. A microscope slide substantially as hereinbefore described with reference to the accompanying drawings.

ABSTRACT

MICROARRAY SLIDES

- 5 An adhesive label 3 for attaching to a microscope slide 1 is disclosed. The label 3 may have one or more grids or arrays 2 provided on it to help an operator to manually spot a microscope slide 1.

[Fig. 1]

